

## 50" Stainless Ductless Fume Hood

Model # SS-350-DCH-ST



The **50" Wide Stainless Ductless Fume Hood** combines powerful air flow with a sturdy and compact enclosure for optimal respiratory and environmental protection. This fume hood helps protect the operator and environment from hazardous fumes or moderately corrosive particulate generated from applications performed within the hood.


The **50" Stainless Ductless Fume Hood's** powerful fan directs fumes into the filter chamber containing high quality filter media selected depending on the application. Filter media choices include **HEPA** [up to 99.97% efficiency on particles 0.3 microns and larger], **ASHRAE** [up to 95% efficiency on particles 0.5 microns and larger], **Activated Carbon**, and **specialty-blended filtration** [*i.e. acid gas, mercury, aldehyde, and ammonia*].

Typical uses for the **50" Stainless Ductless Fume Hood** include *chemical fumes, pharmaceutical compounding, soldering applications, dust removal, biological applications, and a variety of applications requiring fume and dust removal*. The **50" Stainless Ductless Fume Hood** can provide an effective and economical solution for many commercial, medical, laboratory, and industrial applications. For additional visibility, a transparent version is available at extra cost.

*Our line of ductless fume hoods is not intended for sterile compounding applications.*



50" Stainless Ductless Fume Hood

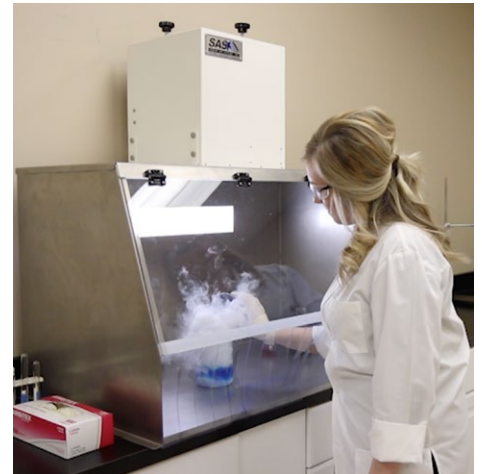

50" DCH provides fume control for various applications.

### Product Features

- **Includes:** variable speed controller (115V only), velometer, and LED light
- **Optional:** clear anti-static vinyl curtains, Magnehelic Gauge, hour counter, and spill tray
- **Standard Sizes:** 30", 40", and 50" (*custom sizes available*)
- Quiet operation
- Long filter life
- Simple, quick "no tool" filter change
- Portable and lightweight
- Reliable, low maintenance operation




21221 FM 529 Rd, Cypress, TX 77433  
1.800.799.4609 • sales@sentryair.com

[www.sentryair.com](http://www.sentryair.com)

# Fume & Particulate Removal

## Product Specifications

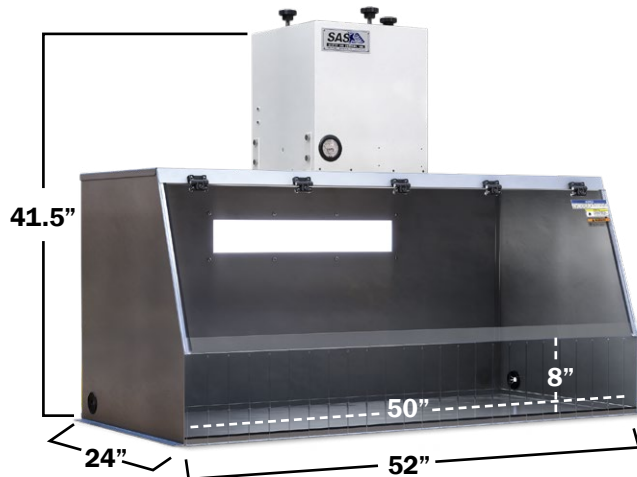

Shown with outer dimensions.  
Dimensions are approximate.  
Shown with optional curtains.



21221 FM 529 Rd  
Cypress, TX 77433  
1.800.799.4609  
sales@sentryair.com

[www.sentryair.com](http://www.sentryair.com)

<b>Hood Dimensions (D x W x H)</b>	<b>Outer:</b> 24" x 52" x 25" <b>Interior:</b> 23.5" x 50" x 24" <b>Access Area:</b> 8" x 50"
<b>Hood Material</b>	<b>Hood:</b> 16 awg 304 stainless steel <b>Fasteners:</b> 304 stainless steel <b>Front Window:</b> 0.25" clear polycarbonate
<b>Blower Dimensions</b>	15.5" L x 13" W x 14.5" H
<b>Blower Material</b>	16 Ga. carbon steel with polyester epoxy coating on all surfaces
<b>Weight</b>	125 lbs
<b>Air Volume (approximately)</b>	350 CFM <i>(varies with filter media)</i>
<b>Inlet Velocities on High (approximately)</b>	<b>Carbon Pre-filter + HEPA Filter</b> – 105 FPM high, 25 FPM low <b>Carbon Pre-filter + Heavy Duty Carbon Filter</b> – 80 FPM high, 20 FPM low
<b>Electrical</b>	115/1/60, 2 amps 220/1/50, 0.75 amps 8' grounded power cord with NEMA 5-15P plug
<b>Sound Level</b>	61 to 66 dba (low to high)
<b>Maximum Temperature</b>	Operating temperature not to exceed 250° F
<b>Main Filter Options</b>	<b>HEPA</b> (up to 99.97% efficient on particles down to 0.3 microns) <b>ASHRAE</b> (up to 95% efficient on particles down to 0.5 microns) <b>Activated Carbon Specialty-Blended Filter Media</b> (i.e. acid gas, mercury, aldehyde, and ammonia)
<b>Warranty</b>	Limited two-year warranty from date of shipment on defects due to materials or workmanship.
<b>Patent #</b>	5,843,197